


Wax Therapy Appliance

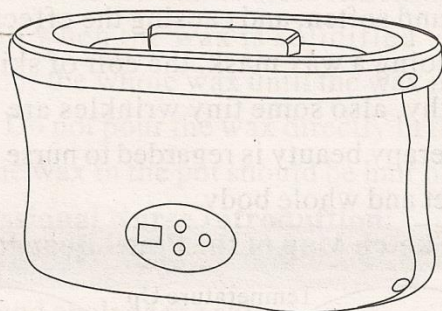


User's Manual

Thanks very much for purchasing our Ion Sprayer. Please read carefully our User's Manual prior to use. After fully mastered it, you will see clearly its excellent performance and operational method.

Part 1. Functions Introduction:

The wax therapy is favorite by beauty shop and families consumers for its remarkable effect and easy to operate. There is good effect for dry nature skin or ageing wrinkles skin, specially in autumn or winter seasons. There are full of rich elements such as vitamin, hydrolysis bone glue albumen, remove wrinkles and brightening in the wax. The wax therapy is called crystal mask. The main function is to stimulate effective medicines get into skin to extend secretion of sweat. That cannot be removed sweat will enter into skin horny substance layer leading it to warm and soften, and causing the effect of skin smooth and tender. After doing a wax mask, the color of skin can change to lubricate and healthy, also some tiny wrinkles are vanished along with it. The wax therapy beauty is regarded to nurse face, meanwhile for head, hands, feet and whole body.

Part 2. Operating Sketch Map of the Panel Board:

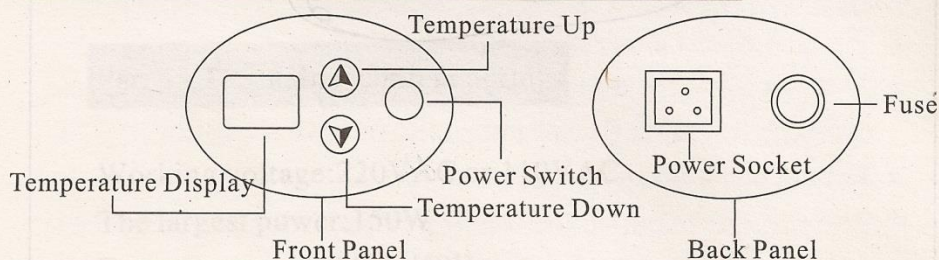

Part 3. Instrument Usage:

1. Plug one end of the power line into the power socket, and connect the other end to commercial current according to the power identification (220V AC / 110V AC).
2. Press down the power switch, then the power is on (with factory setting being 50°C) and the machine starts to work.

3. Adjust the temperature to 65°C for heating ; after the beautifying wax is wholly melted, adjust the temperature to 55°C; the instrument may be used after the wax temperature drops to 45-55°C and keeps within this range.
4. The instrument has the function of keeping constant temperature, and remains people's temperature setting.

Part 4. Points for Attention:

1. Beautician must brush a little wax on the wrist side to make sure the temperature of wax is suitable, avoiding scalding skin before rub the wax liquid on the skin surface.
2. Do not use when skin is damaged.
3. Do not use when there is acne on the skin surface.
4. Please shut down the mains after finishing. Put through the mains again when the wax is solidified completely. Then does not pour the whole wax until the wax near the pot surface is melted. Do not pour the wax directly in the liquid situation.
5. Depth of the wax in the pot should be half height of the pot.

Part 5. Professional Nurse Introduction:

Professional face nurse

1. Wash face and analyze the skin.
2. Steaming face for the minutes rub massage lotion on the face, and massage for about five minutes, then remove the lotion.
3. Rub on the tending lotion.
4. Rub the wax mask on the tending lotion if the temperature is suitable after testing in wrist side.

Professional hands nurse

1. Wash hands and remove the ornament in the hands. Test the max temperature if suitable with your finger.
2. Spread the fingers, soak in the wax liquid for 2-3 seconds then take them out.

3. Duplicate the following operation for 4-6 times. At this moment, hands are solidified to thick wax mask.
4. Wrap up the hands with gluey gloves right now. So that can maintain the temperature and absorb nourishment for pore.
5. After 15-20 minutes, remove the gluey gloves, then remove the wax mask and extra fat of the fingers.

Hands will be come to delicate and smooth.

Professional face nurse

1. Wash feet and test the wax temperature if suitable with your finger.
2. Soak the feet in the wax liquid for 2-3 seconds then take out (Rub the wax on the feet with brush if the feet can't soak in the pot.)
3. Duplicate the following operation for 4-6 times. At this moment. Feet are solidified to thick wax mask.
4. Wrap up the feet with gluey bag right now so that can maintain the temperature and skin can be heated completely.
5. After 15-20 minutes. Remove the gluey bag. Then remove the wax mask and extra fat of the foot-nail, feet can be come to delicate and springy.

Part 6. The main technical norm:

Working voltage: 220VAC or 110VAC

The largest power: 150W

Frequency range: 50Hz/60Hz

Attention:

The appliance has been adjusted and inspected strictly before offering for sale. Must repaired by the agency and special center in every place of our corporation if breakdown during use. Do not open the machine body to repair by unprofessional staff.